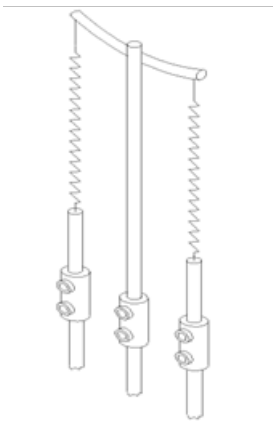

Filament Replacement

Tungsten filaments have an emission life of one to three years in UHV. The life will be dependent on the residual gases present, range of pressures and vibration coupled to the filament. Tungsten filaments become brittle in use and the most common causes of failure are mechanical or thermal shocks. If the gauge head is to be moved to another location it is advisable to replace the filament assembly.

Replacement filament assemblies are available from AML distributors at low cost. The filament assemblies are fixed to the feedthrough assembly by barrel connectors and Allen-headed set screws. An Allen key is provided with each replacement filament assembly. A pair of long-nosed pliers, side cutters, cleanroom gloves, cotton buds and some isopropyl alcohol are required.



Clean the jaws of the cutters and pliers with alcohol. Swab some alcohol onto the set screws to lubricate them. Support the feedthrough wire with the long-nosed pliers close to the barrel connector and loosen the screw holding the filament support. This is the screw furthest away from the feedthrough flange. Loosen the screws on all three supports and remove the old filament assembly. Remove the restraining sleeve from the new assembly. Fit the centre support of the new assembly into the barrel connector and tighten the screws. Hold each of the filament terminations in turn with the pliers and extend them down to the outer barrel connectors. Fix them in place with the set screws. Ensure that the centre support of the filament assembly is parallel to the gauge axis and approximately 3mm from the side of the grid.
